


UKMED/CAL®

*It's Interventional.*

# *It's the Neo-Flex Ureteric Stent/*

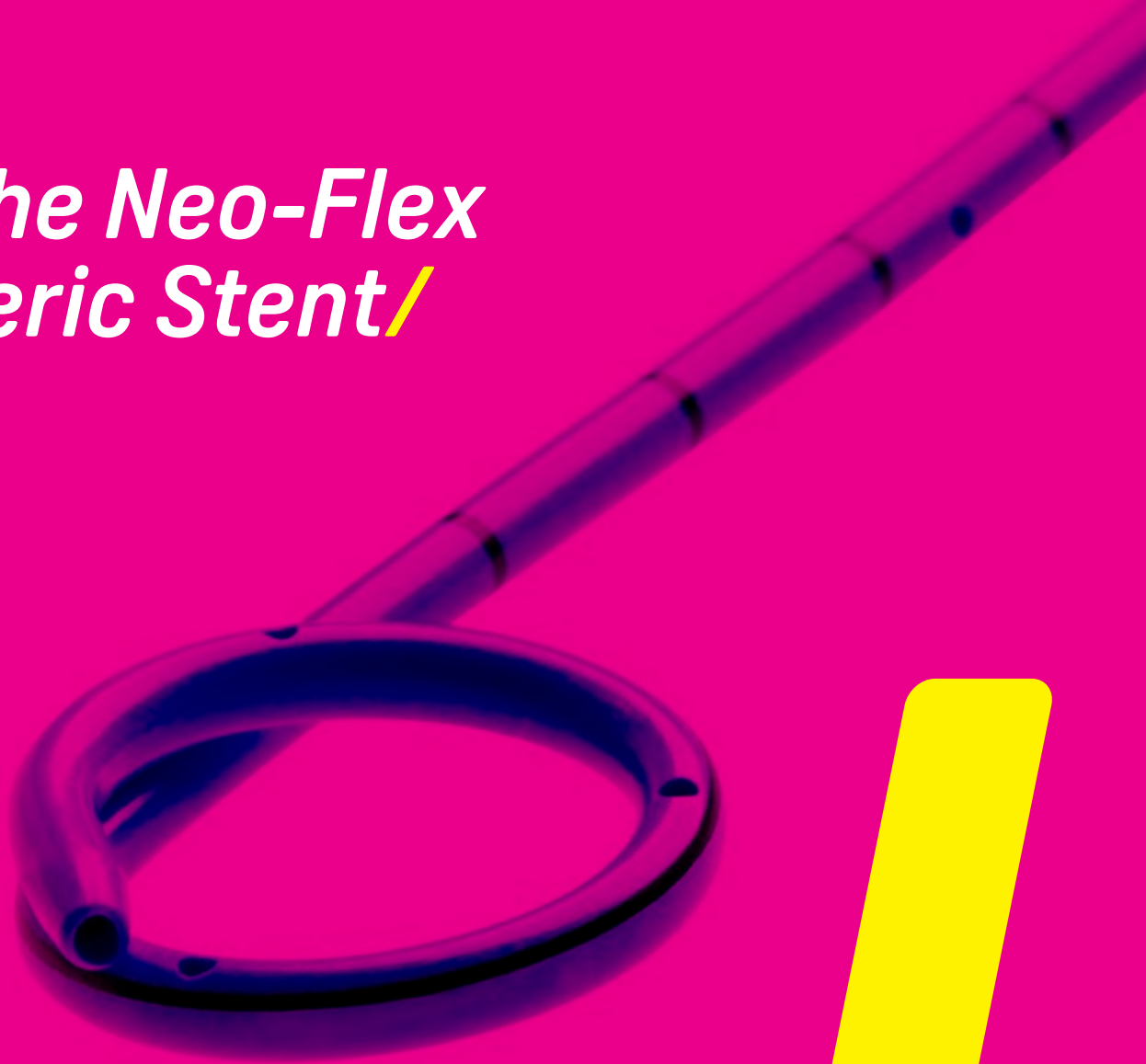

@UK\_Medical  
hello@ukmedical.com  
ukmedical.com



# *It's about clinical flexibility... /*

*The Neo-Flex ureteral stent's enhanced durometer material provides strong resistance to crush related occlusions, whilst simultaneously providing a high degree of patient comfort.*


*Perfectly suited for both short and long-term placement, Neo-Flex is available in a wide range of sizes and is designed to provide all round clinical flexibility.*

## **Product benefits:**

- ***Provides smooth, controlled advancement through the renal structure***
- ***Suitable for antegrade and retrograde procedures***
- ***Precise fluoroscopic visualisation with reassured stent position***
- ***Reliable, long-term stent patency***
- ***Ideally suited for malignant strictures, normally associated with rapid stent occlusion***

## **Additional information?**

If you would like further assistance with using our ureteric stents, or need advice on how to order, please call us on: **0114 268 8880**

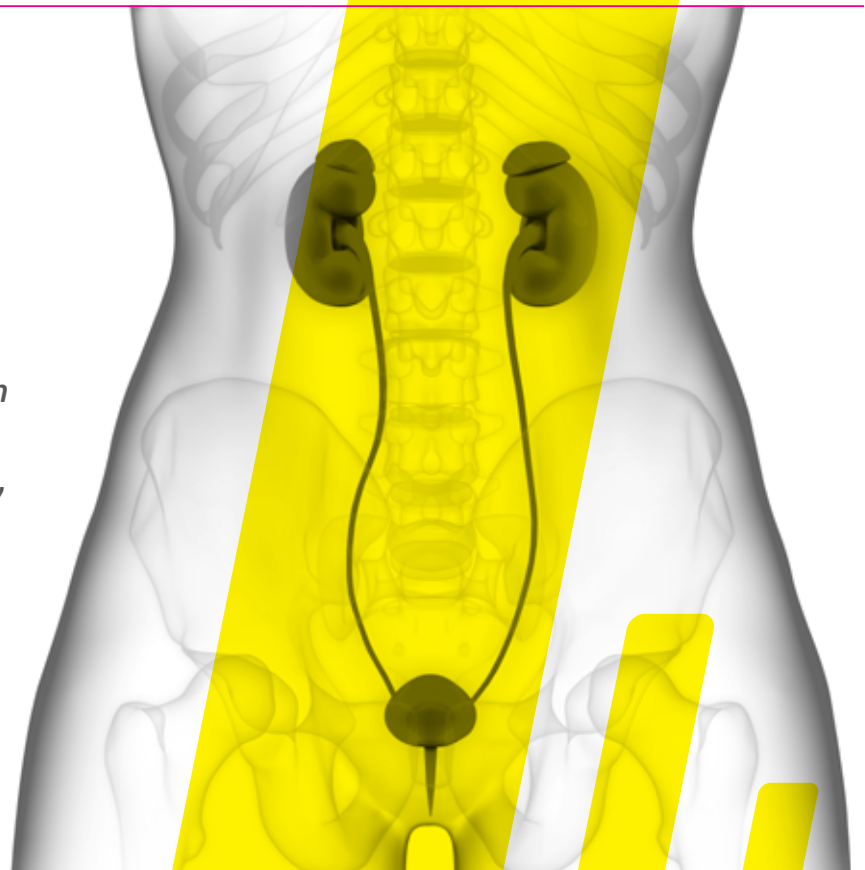

**Product features:**


Enhanced durometer material

High quality hydrophilic coating

Smooth tapered tip


Optimal balance between barium sulphate content and flexibility

The thin-walled Neo-Flex has a large lumen, providing an optimal ID:OD ratio

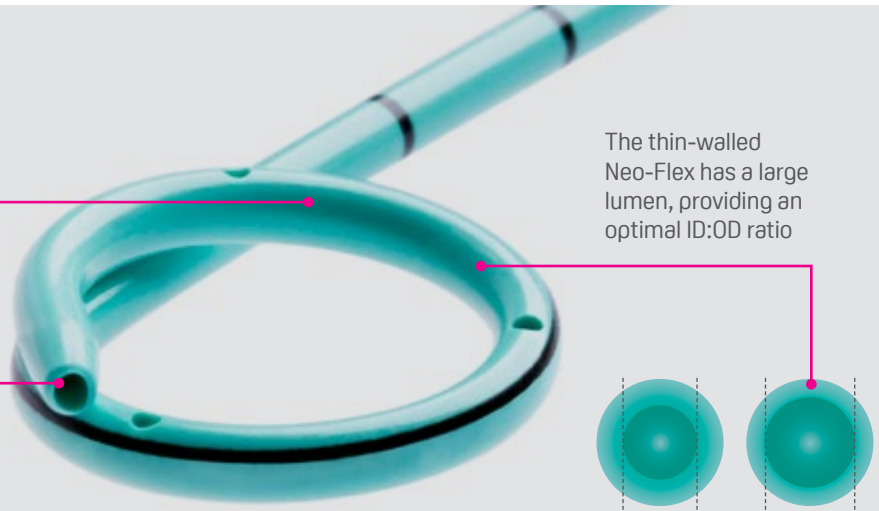

Skived drainage holes ensuring a smooth surface, free of luminal artefact and/or porosity

Six offset drainage holes around the pigtails with spiralled drainage holes along the working length

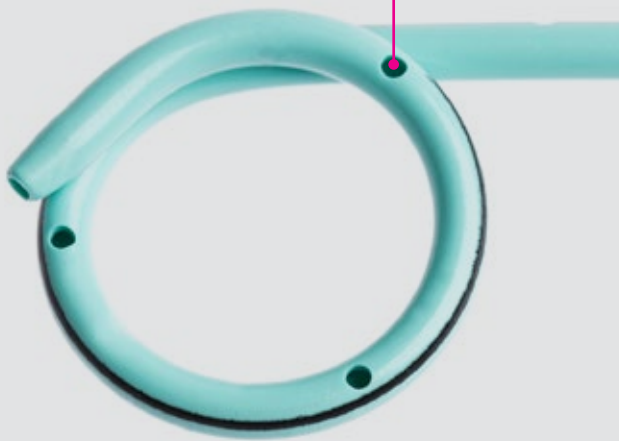

Funnelled profile of skived drainage holes helps prevent encrustation and maintain reliable flow rates

Eliminates artefact for a smoother internal edge

Smoother hole surface/cross section

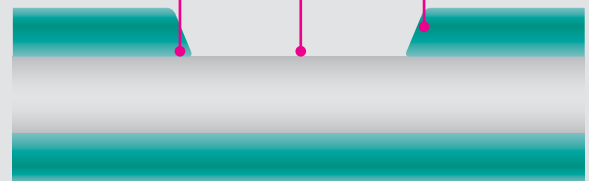


Fig 1. Cross section illustration of skived drainage holes

Smooth coated pusher/positioner with radiopaque marker



Rapid pigtail formation in the renal pelvis and bladder following guidewire removal

Centimetre markings along the working length of the stent



## Kit contents:

- Pusher with radiopaque marker
- Clamp
- 0.035 guidewire (fixed core PTFE coated)
- Double J stent with retrieval thread



## Stocked Items

Product Code	Description	Boxed
BTDJ6024-211-HSL	Ureteric stent 6Fr x 24cm	5
BTDJ6028-211-HSL	Ureteric stent 6Fr x 28cm	5
BTDJ7024-211-HSL	Ureteric stent 7Fr x 24cm	5
BTDJ7028-211-HSL	Ureteric stent 7Fr x 28cm	5
BTDJ8024-211-HSL	Ureteric stent 8Fr x 24cm	5
BTDJ8028-211-HSL	Ureteric stent 8Fr x 28cm	5
BTDJ8030-211-HSL	Ureteric stent 8Fr x 30cm	5

Also available in the following sizes:  
4 to 8 French and lengths 20 to 30cm  
in two-centimetre increments.

### Also available:

- Pyelo open/closed
- Mono J
- Multi-length
- Dual durometer

Manufactured by

**BIOTEQ**<sup>®</sup>

Distributed by:

UK Medical  
Albreda House  
Lydgate Lane  
Sheffield, S10 5FH

T: +44 (0) 114 268 8880  
F: +44 (0) 114 268 8881  
sales@ukmedical.com  
www.ukmedical.com

**UKMED/CAL**<sup>®</sup>  
It's Interventional.