


UKMED/CAL®

It's Interventional.

It's the Nitinella Plus Biliary Stent/

hello@ukmedical.com
ukmedical.com

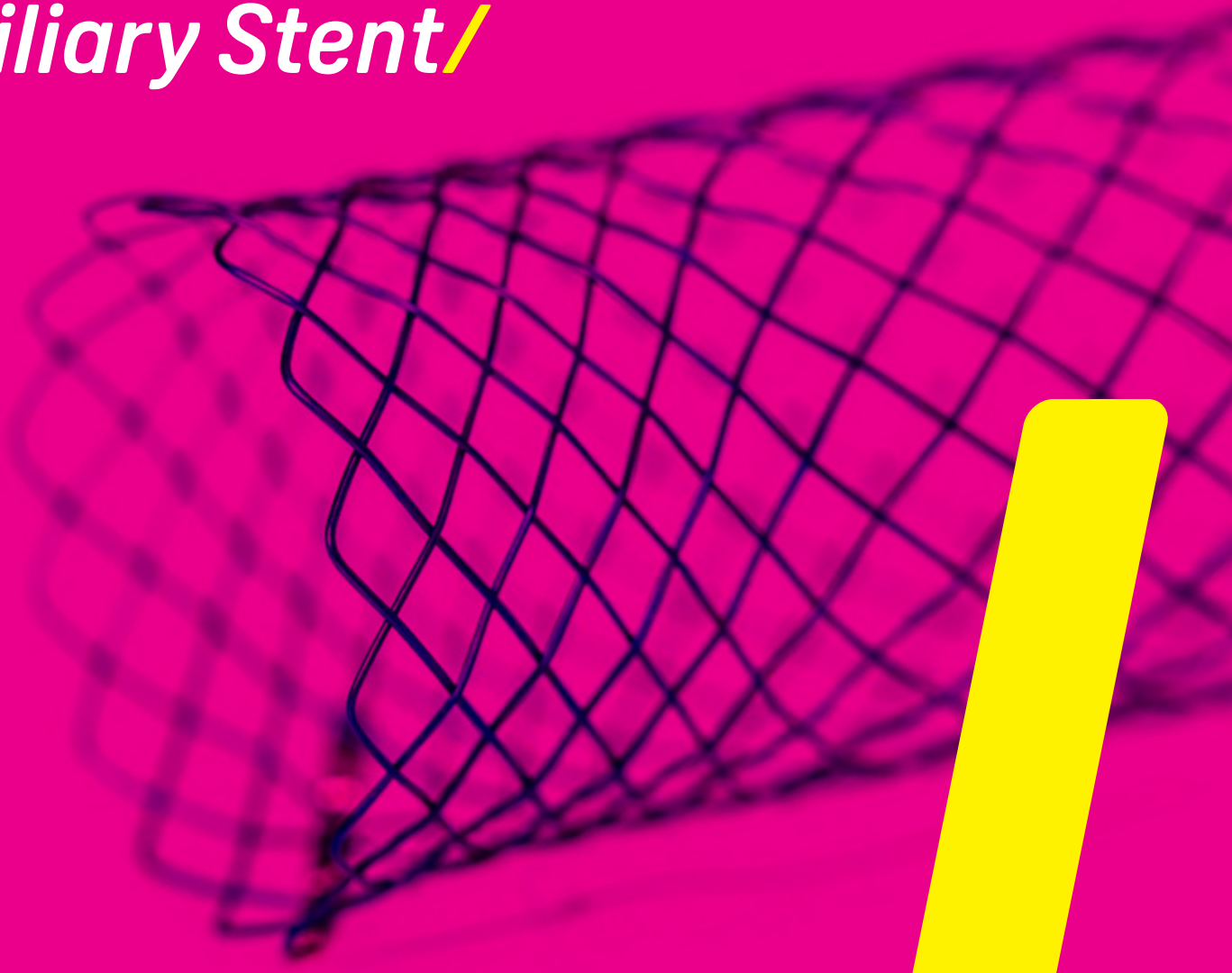

It's in the design.../

The Nitinella Plus Biliary Stent provides an optimal balance between patency, conformability and a reduced risk of migration.

The atraumatic stent is suitable for both transhepatic and endoscopic implant techniques.

Available in uncovered, partially covered or fully covered versions, Nitinella Plus provides a variety of options to address the demands posed by malignant strictures whilst helping to improve patient outcomes.

Product benefits:

- ***Precise fluoroscopic visualisation with reassured stent position***
- ***Reliable long-term stent patency and patient comfort***
- ***Minimised risk of migration reduces requirement for repeat procedures***
- ***Suitable for transhepatic and endoscopic implant techniques***
- ***Reduced stent shortening upon release aids accurate placement***
- ***Minimised risk of perforation***

Additional information?



If you would like further assistance with using our Nitinella Plus biliary stent, or need advice on how to order, please call us on: **0114 268 8880**




Product features:

Variable weave at distal and proximal ends provides additional radial force to minimise stent migration


Enhanced radial force maintains stent patency

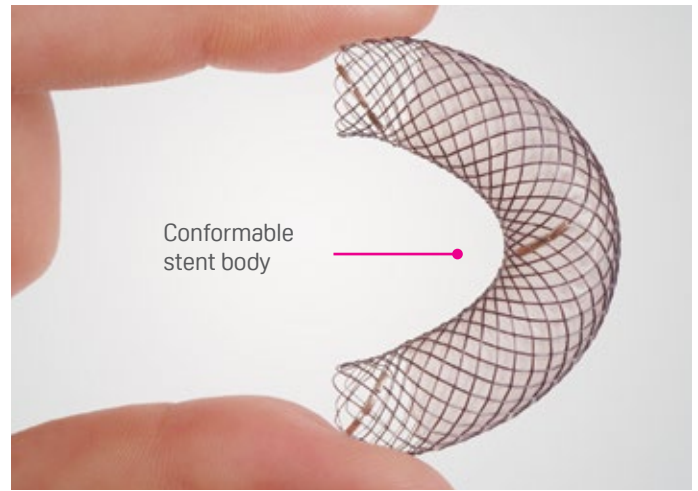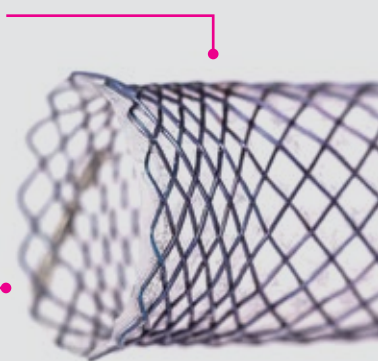


Conformable stent body

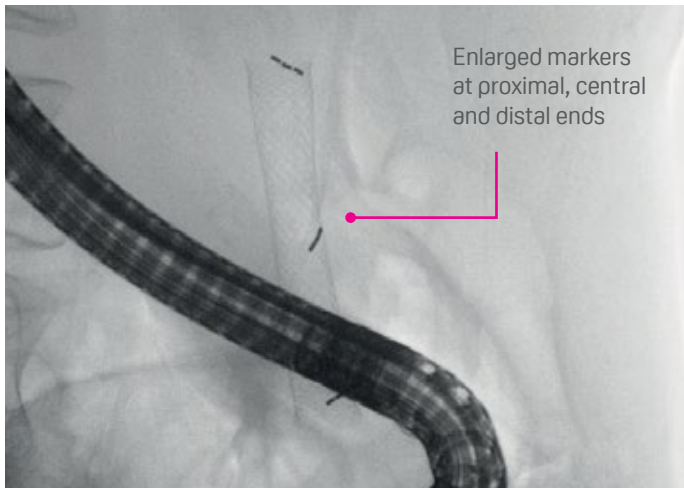


Enlarged markers at proximal, central and distal ends

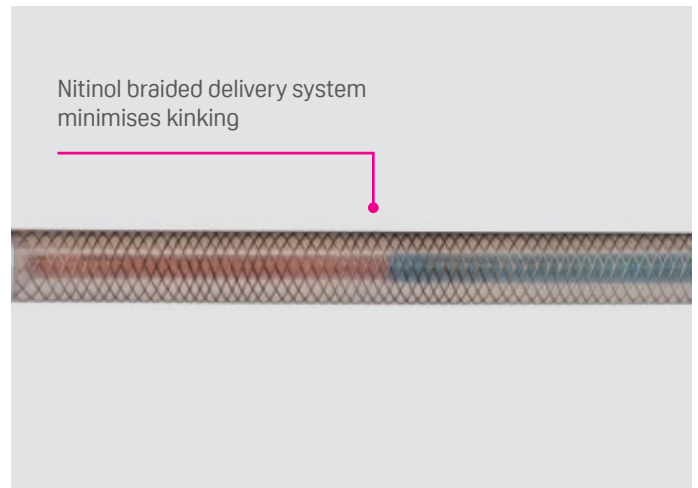


Nitinol braided delivery system minimises kinking

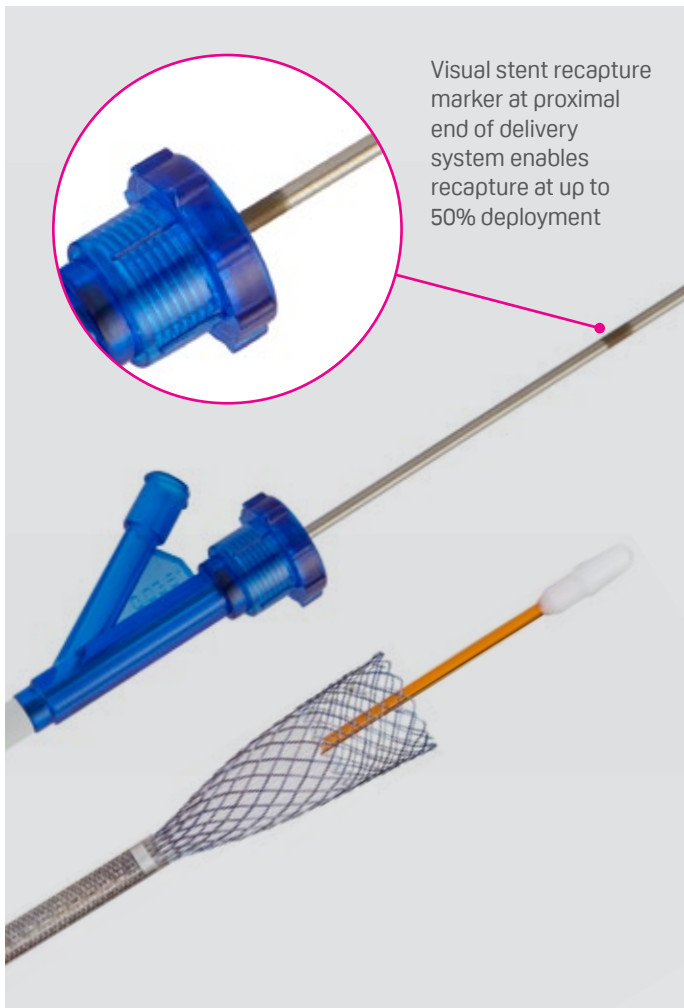

Visual stent recapture marker at proximal end of delivery system enables recapture at up to 50% deployment

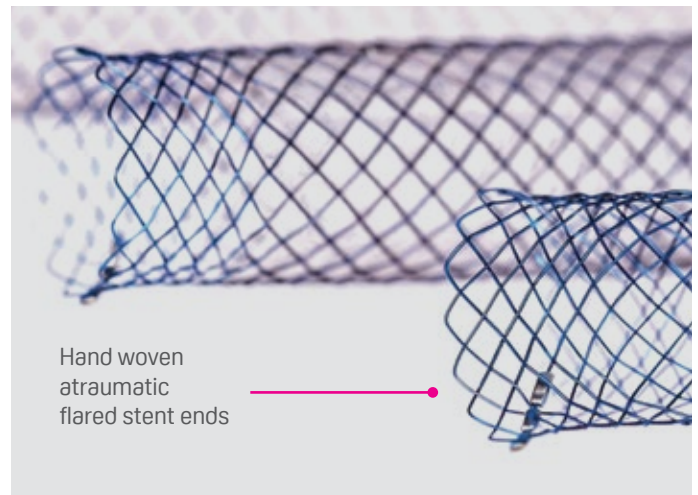

Hand woven atraumatic flared stent ends

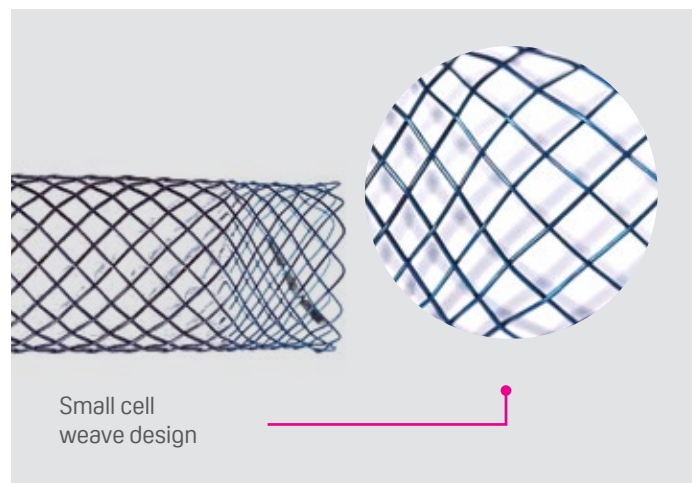

Small cell weave design

Transhepatic, 60cm delivery system

Stent Diameter (mm)	Nominal Length (mm)	Delivery System diameter (F)	Product Code Uncovered	Product Code Partially Covered	Product Code Fully Covered
8	40	7 for uncovered, 9 for partially & fully covered	019-01B-08-040-T	019-01C-08-040-T	019-01D-08-040-T
	60		019-01B-08-060-T	019-01C-08-060-T	019-01D-08-060-T
	80		019-01B-08-080-T	019-01C-08-080-T	019-01D-08-080-T
	100		019-01B-08-100-T	019-01C-08-100-T	019-01D-08-100-T
10	40	7 for uncovered, 9 for partially & fully covered	019-01B-10-040-T	019-01C-10-040-T	019-01D-10-040-T
	60		019-01B-10-060-T	019-01C-10-060-T	019-01D-10-060-T
	80		019-01B-10-080-T	019-01C-10-080-T	019-01D-10-080-T
	100		019-01B-10-100-T	019-01C-10-100-T	019-01D-10-100-T

0.035" x 180cm ultra-stiff guidewire recommended

Endoscopic, 175cm delivery system

Stent Diameter (mm)	Nominal Length (mm)	Delivery System diameter (F)	Product Code Uncovered	Product Code Partially Covered	Product Code Fully Covered
8	40	7 for uncovered, 9 for partially & fully covered	019-01B-08-040	019-01C-08-040	019-01D-08-040
	60		019-01B-08-060	019-01C-08-060	019-01D-08-060
	80		019-01B-08-080	019-01C-08-080	019-01D-08-080
	100		019-01B-08-100	019-01C-08-100	019-01D-08-100
10	40	7 for uncovered, 9 for partially & fully covered	019-01B-10-040	019-01C-10-040	019-01D-10-040
	60		019-01B-10-060	019-01C-10-060	019-01D-10-060
	80		019-01B-10-080	019-01C-10-080	019-01D-10-080
	100		019-01B-10-100	019-01C-10-100	019-01D-10-100

0.035" x 360cm ultra-stiff guidewire recommended

Related Products

Product Code	Description	Box Qty
BT-PDS-XXXX-W-RB*	Biliary catheter int/ext	5
BT-PDN-618-RB	Percutaneous introducer set w/o GW	5
BT-PDN-618-RB+GW2	Percutaneous introducer set 90cm GW	5
BT-PDN-618-RB+GW2A	Percutaneous introducer set 150cm GW	5

*XXXX denotes length and French size, please specify when ordering.

Full range of GI stents available. Please contact us for more information.

Manufactured by



Distributed by:

UK Medical
Albreda House
Lydgate Lane
Sheffield, S10 5FH

T: +44 (0) 114 268 8880
F: +44 (0) 114 268 8881
sales@ukmedical.com
www.ukmedical.com

@UK_Medical

UKMEDICAL®
It's Interventional.