

**UKMED/CAL®**

*It's Interventional.*

# *It's Nephrostomy Solutions/*

*hello@ukmedical.com*  
*ukmedical.com*



# *It's about the complete solution... /*

*UK Medical's aim is to provide a complete solution for your nephrostomy procedures.*


*Designed with patient comfort and procedural ease in mind our comprehensive nephrostomy package helps to simplify your ordering process and streamline the patient care pathway.*

## **Product benefits:**

- ***Complete nephrostomy procedure solution***
- ***Catheter is suitable for long-term placement***
- ***Provides atraumatic, controlled advancement through the renal structure***
- ***Designed for enhanced patient comfort***
- ***Secure and comfortable catheter fixation***
- ***Hospital to home solution***

## **Additional information?**

If you would like further information about Nephrostomy Solutions, or need advice on how to order, please call us on: **0114 268 8880**

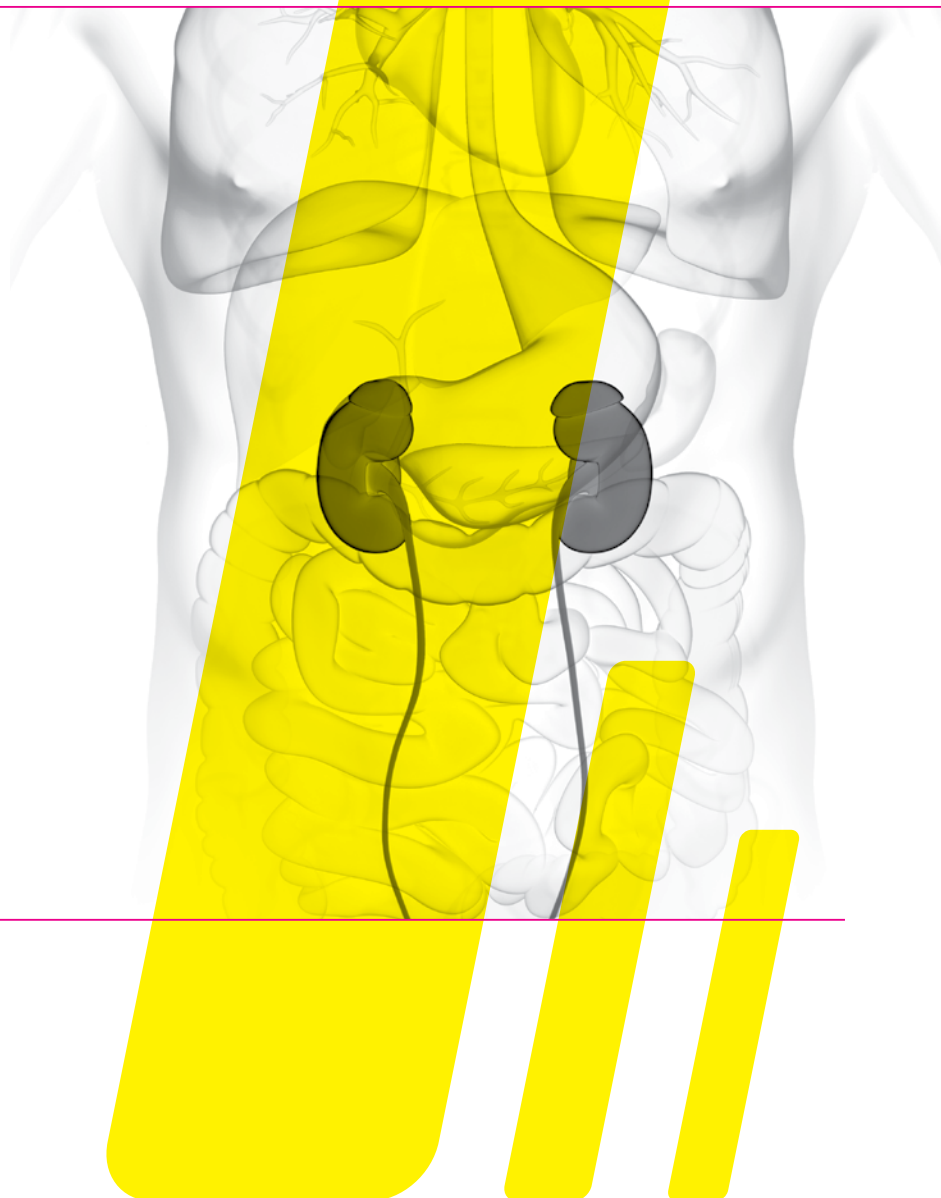


**Product features:**

Smooth tapered plinth tip

Large skived drainage holes, highly resistant to encrustation

Multiple pigtail design options available on request

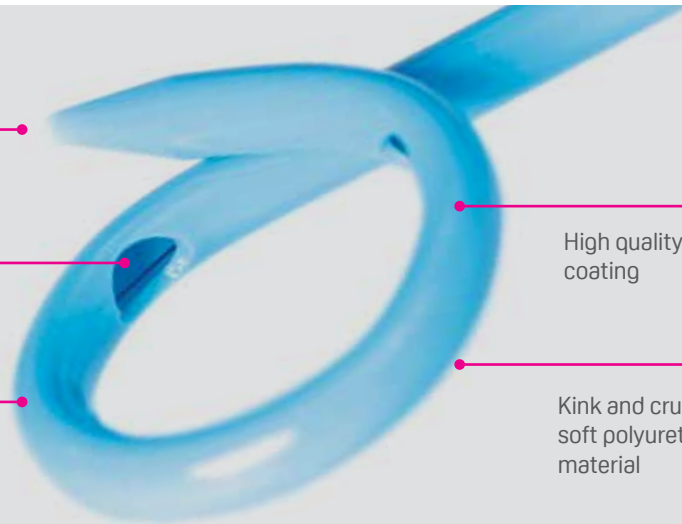

High quality hydrophilic coating

Kink and crush resistant soft polyurethane material

**Choice of catheter fixation devices available:**

**Drain Guard**



- Transparent hydrocolloid cover allows visualisation of puncture site
- Hydrocolloid adhesive
- Catheter support block designed to prevent kinking

**Drain Lok**



- Alternative catheter fixation to suit your patients' needs
- Hydrocolloid adhesive
- Catheter support block designed to prevent kinking

**Introducer**




Echogenic 21G introducer needle with bevel tip



3 part dilator with radiopaque marker



0.018 nitinol fine guidewire with palladium tip

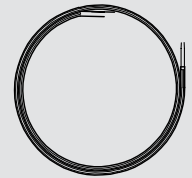

Optional 0.035 guidewire

**Drainage Bags**

500ml to 2L options available



## Ordering Information

### Nephrostomy Drainage Catheters

Description	Design type	Product code	Size	Working length	Box quantity
Neo-Hydro Drainage Catheter - LOCKING PIGTAIL	Standard Pigtail	BT-PD1-0630-W	6Fr x 30cm	29cm	5
	Standard Pigtail	BT-PD1-0640-W	6Fr x 40cm	39cm	5
	Standard Pigtail	BT-PD1-0730-W	7Fr x 30cm	29cm	5
	Standard Pigtail	BT-PD1-0820-W	8Fr x 20cm	19cm	5
	Standard Pigtail	BT-PD1-0825-W	8Fr x 25cm	24cm	5
	Standard Pigtail	BT-PD1-0830-W	8Fr x 30cm	29cm	5
	Standard Pigtail	BT-PD1-0840-W	8Fr x 40cm	39cm	5
	Standard Pigtail	BT-PD1-1030-W	10Fr x 30cm	29cm	5
	Standard Pigtail	BT-PD1-1230-W	12Fr x 30cm	29cm	5
	Standard Pigtail	BT-PD1-1430-W	14Fr x 30cm	29cm	5
Neo-Sof Dedicated Nephrostomy Drainage Catheter - LOCKING PIGTAIL	Inc. Flexible Stiffening Cannula	BT-PD1-0830-W-FSC	8Fr x 30cm	29cm	5
	Mini-Closed Loop/FSC/No Stylet	BT-PDS-0830-W-MN	8Fr x 30cm	29cm	5
	Standard Closed Loop	BT-PDS-0830-W-CL	8Fr x 30cm	29cm	5

Additional sizes, pigtail designs and catheter compositions available on request.

### Accessory Products

Description	Design type	Product code	Box quantity
Percutaneous introducer set	Standard Introducer kit	BTPDN618RB	5
Percutaneous introducer set 90cm GW.	Std kit Inc. 90cm 035 GW	BTPDN618RB+GW2	5
Percutaneous introducer set 150cm GW.	Std kit Inc. 150cm 035 GW	BTPDN618RB+GW2A	5
Drainage bag 500ml	Inc. ML Luer Lock Drainage Line	LU500	25
Drainage bag 800ml	Inc. ML Luer Lock Drainage Line	LU800	25
Drainage bag 2L	Inc. ML Luer Lock Drainage Line	UB-121T-1(ML)	25
Leg bag with fabric backing and velcro strap 800ml	Inc. ML Luer Lock Drainage Line	DB-800C-MS-V	10
Drain-Guard fixation device	Oval - Anti-Kink - Hydrocolloid	AD-8150	25
Drain-Lok fixation device	Square - Anti-kink - Hydrocolloid	BT-RD-DL	25

Further guidewire and drainage bag options available on request.

Distributed by:

UK Medical  
Albreda House  
Lydgate Lane  
Sheffield, S10 5FH

T: +44 (0) 114 268 8880  
F: +44 (0) 114 268 8881  
sales@ukmedical.com  
www.ukmedical.com

 @UK\_Medical

 **UKMEDICAL**<sup>®</sup>  
It's Interventional.