


UKMEDICAL®

It's Interventional.

*It's the
Danis stent/*

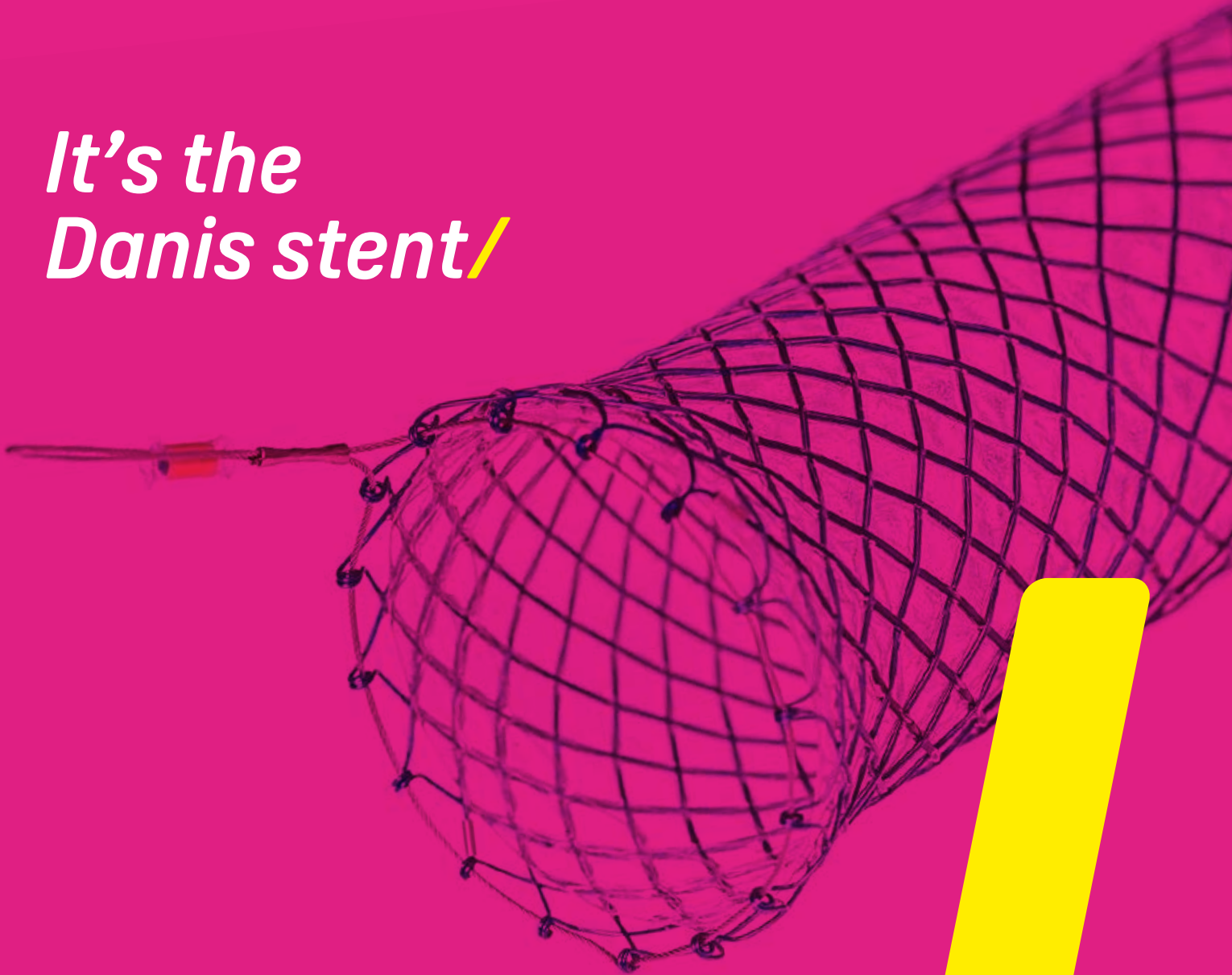

*For oesophageal
varices/*

*hello@ukmedical.com
ukmedical.com*



It's all about patient outcomes... /

The Ella Danis stent provides a fast, safe and effective treatment option for patients with refractory oesophageal variceal bleeds.


Designed to tamponade oesophageal varices whilst maintaining normal dietary intake and minimising the need for high-dependency observation, Danis can stabilise your patient allowing time to plan future treatment options and improve patient outcomes.

Product benefits:


- **Improved patient quality of life compared to balloon tamponade**
- **A 2010 study concluded, 'insertion of the SX-Ella Danis stent in patients with refractory variceal bleeding or complications of previous therapy is effective for the control of bleeding'¹**
- **Reduces potential resource-drain on high dependency units**
- **Can be placed without imaging guidance in emergency situations**
- **Significantly reduced complications and SAEs when compared to balloon tamponade in acute oesophageal varices²**

Additional information?

If you would like further assistance with using our Danis stent, or need advice on how to order, please call us on: **0114 268 8880**

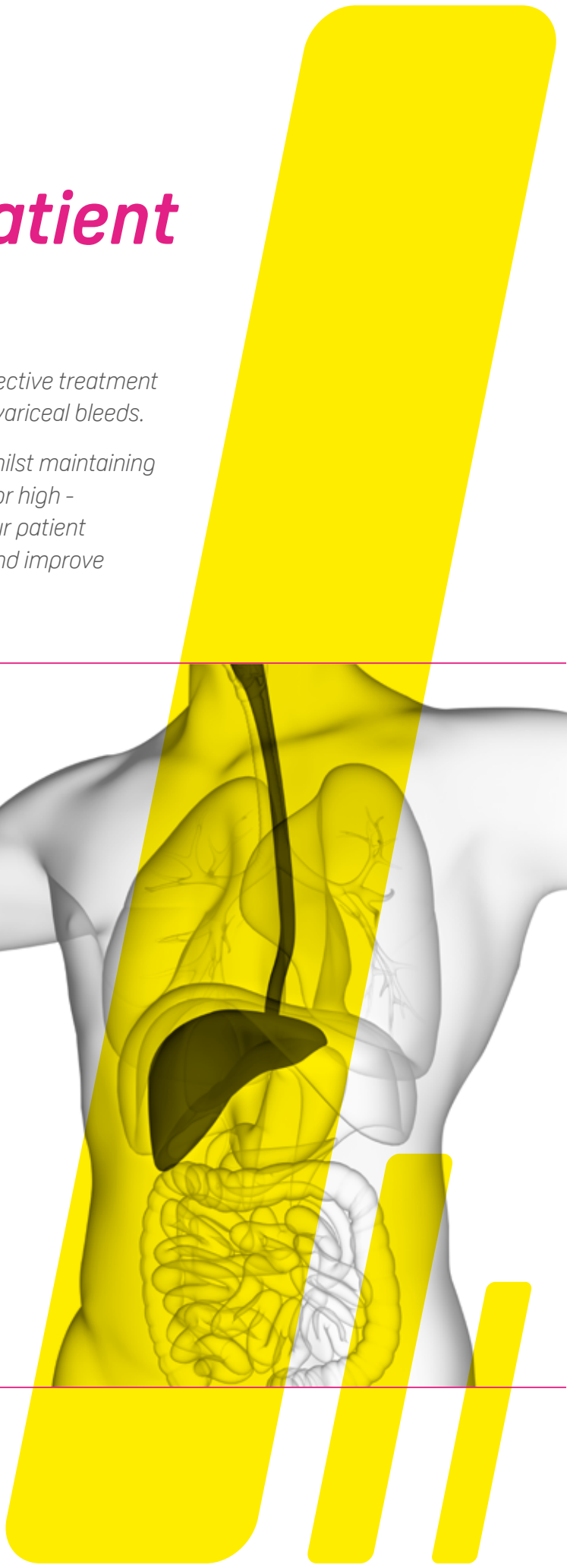


Product features:

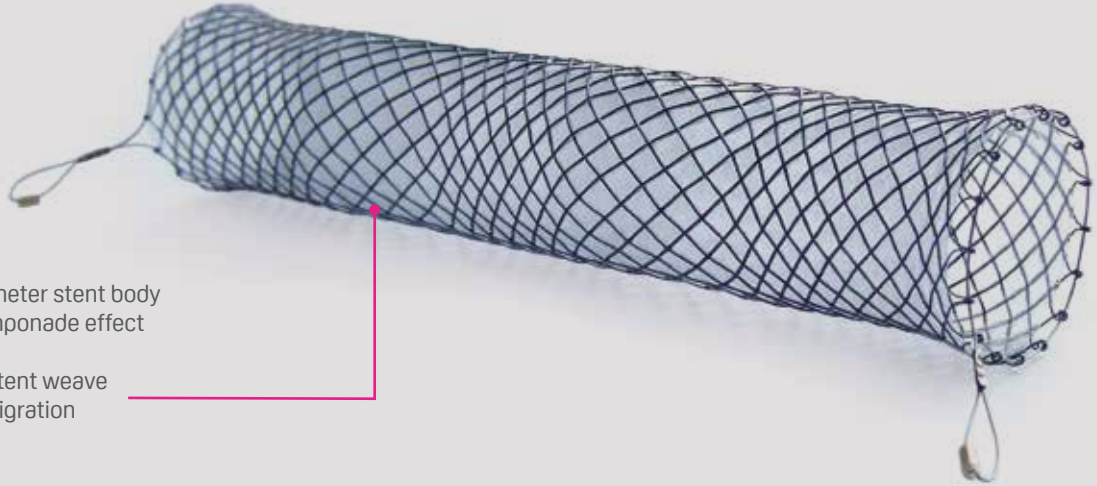


Large 25mm diameter stent body for enhanced tamponade effect

Unique variable stent weave reduces risk of migration

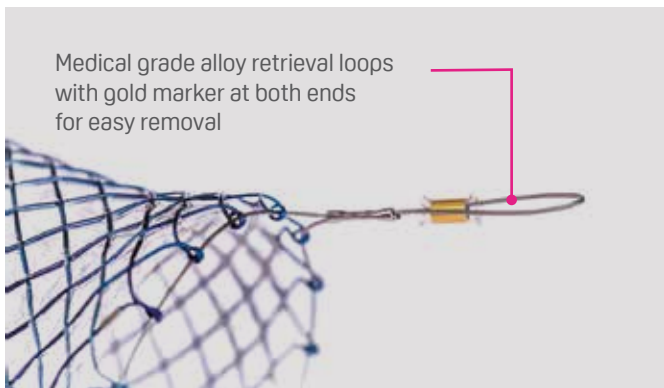


Medical grade alloy retrieval loops with gold marker at both ends for easy removal

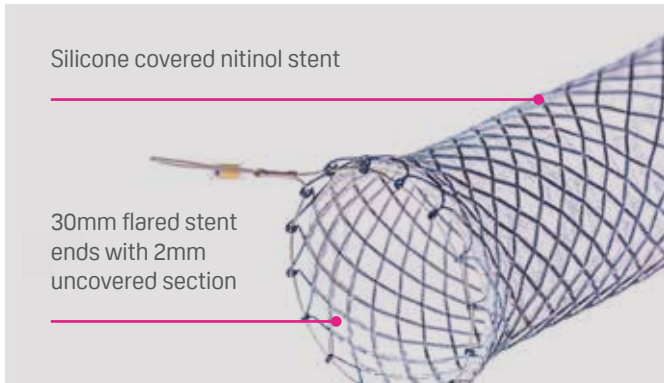

Silicone covered nitinol stent

30mm flared stent ends with 2mm uncovered section

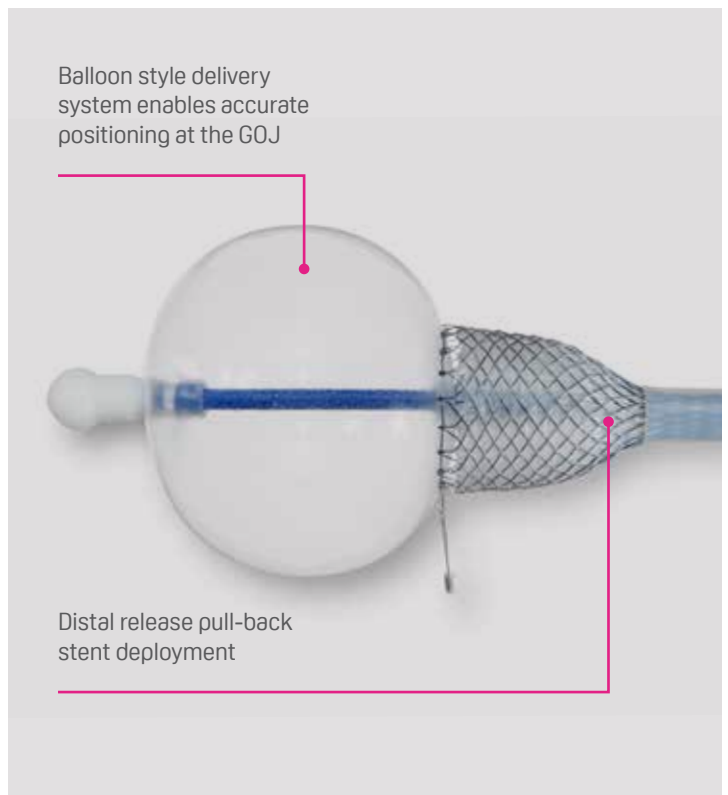


Balloon style delivery system enables accurate positioning at the GOJ

Distal release pull-back stent deployment

It's the Results²


Control of bleeding

- Danis
- Sengstaken-blakemore

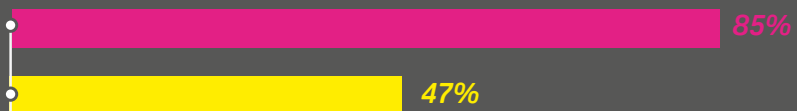

Serious Adverse Events

- Danis
- Sengstaken-blakemore

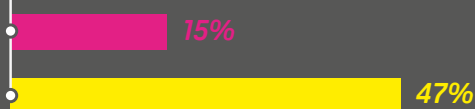

Transfusion

- Danis patients
- Sengstaken-blakemore patients



Kit contents:

- 28F delivery system with pre-loaded stent
- Balloon valve
- 50ml syringe
- 0.035 guidewire




Ordering Information

Product	Product Code	Nominal Diameter (mm)	Nominal Length
Danis Stent	019-08S-25-135-B	25	135mm
ELLA Extractor System	015-01-28-1000	9.4	100cm

Additional information:

- Delivery system 28F / 20F, active length 60cm
- Ella extractor assists with safe, atraumatic removal of the Danis stent

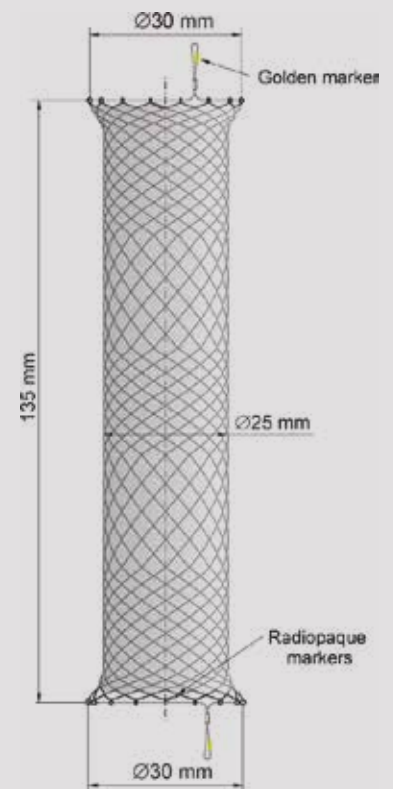

Manufactured by



References:

1. Wright GI, Lewis H, Hogan B, Burroughs, A, Patch D, O'Beirne J, 'A self-expanding metal stent for complicated variceal hemorrhage' (2010) Gastrointestinal Endoscopy
2. Escorsell A, et al., 'Esophageal Balloon Tamponade Versus Esophageal Stent in Controlling Acute Refractory Variceal Bleeding: A Multicenter Randomized Controlled Trial' (2016) Hepatology, (Vol 63, No.6)

Distributed by:

UK Medical
Albreda House
Lydgate Lane
Sheffield, S10 5FH

T: +44 (0) 114 268 8880
F: +44 (0) 114 268 8881
sales@ukmedical.com
www.ukmedical.com

@UK_Medical

UKMEDICAL®
It's Interventional.