

UKMED/CAL®

It's Interventional.

It's the Neo-Sof Biliary Catheter/



*@UK_Medical
hello@ukmedical.com
ukmedical.com*



It's about precise positioning... /

The Neo-Sof biliary catheter is designed for ease-of-use, precise positioning and patient comfort. It is available in locking or non-locking pigtail and comes with both a metal and plastic stiffening cannula for procedural flexibility.


The catheter can be supplied with a direct-connection drainage bag that maintains a safe closed system whilst eliminating the need for, and additional cost of an extra connector.

Product benefits:

- ***Excellent visualisation during placement with fluoroscopic guidance***
- ***Facilitates fast drainage***
- ***Minimal migration facilitated by design ensures best patient comfort and positive outcomes***
- ***Atraumatic catheter introduction***
- ***Easy catheter advancement during insertion***
- ***Fewer interventions required due to high resistance to crush related occlusions***
- ***Resistant to encrustation facilitating easy removal***

Additional information?

If you would like further assistance with using our biliary catheters, or need advice on how to order, please call us on: **0114 268 8880**

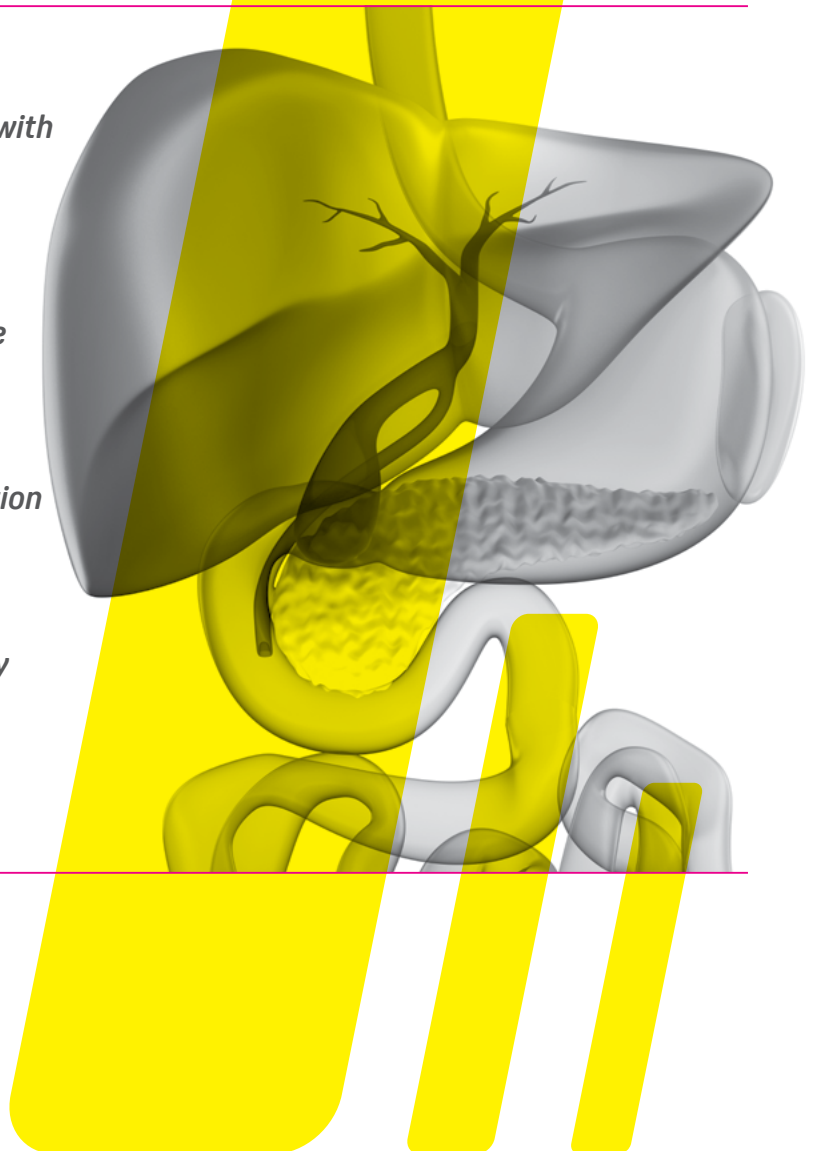


Product features:

Optimal barium sulphate content

High quality hydrophilic coating

Large shaft mounted drainage holes for fast internal & external drainage

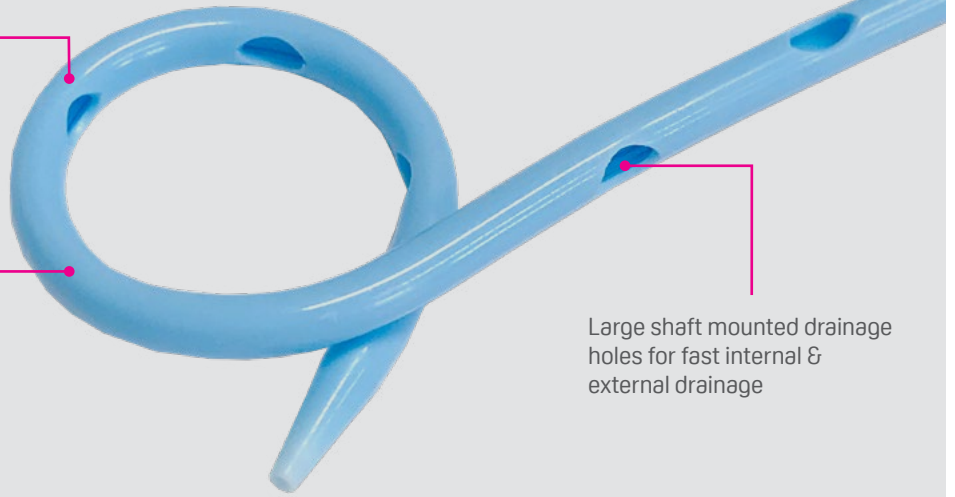

Skived drainage holes with funnelled profile for reduced encrustation



Fig 1. Cross section illustration of skived drainage holes


Enhanced echogenic surface on metal cannula



Splitable pigtail straightener




Radiopaque band indicates final proximal side port

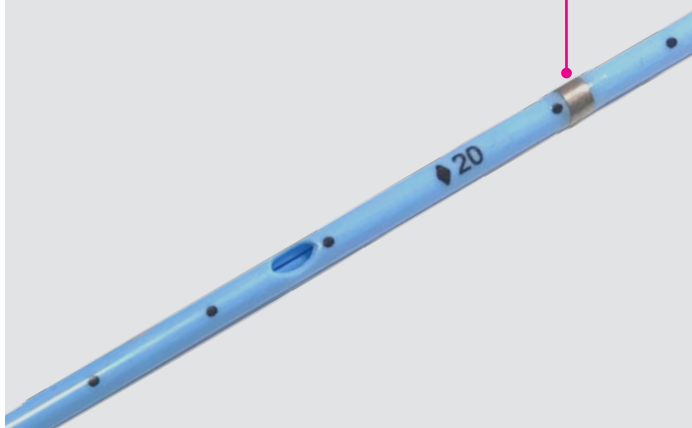

Anti-leak hub design on locking catheter

Catheter size clearly visible on hub

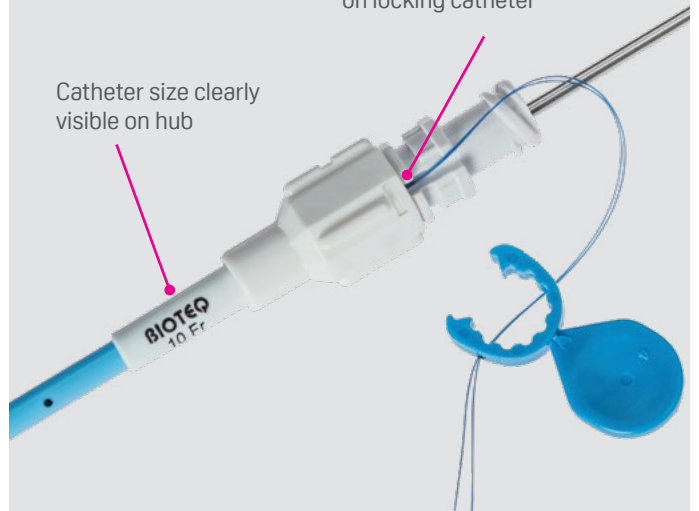

Kit contents:

- Metal stiffener/cannula
- Flexible plastic cannula
- Neo-Sof Biliary Catheter
- Splitable pigtail straightener



Ordering Information

Stocked Items

Description	Diameter	Product code	Length	Box quantity
Neo-Sof Biliary Catheter - Locking	6F	BT-PDS-0640-W-BRB	40cm	5
		BT-PDS-0830-W-BRB	30cm	5
	8F	BT-PDS-0840-W-BRB	40cm	5
		BT-PDS-0850-W-BRB	50cm	5
	10F	BT-PDS-1040-W-BRB	40cm	5
		BT-PDS-1050-W-BRB	50cm	5
	12F	BT-PDS-1240-W-BRB	40cm	5
		BT-PDS-1250-W-BRB	50cm	5

Also available in the following sizes:
6 to 14 French and lengths 20 to 50cm
in five-centimetre increments.

Manufactured by

BIOTEQ[®]

Distributed by:

UK Medical
Albreda House
Lydgate Lane
Sheffield, S10 5FH

T: +44 (0) 114 268 8880
F: +44 (0) 114 268 8881
sales@ukmedical.com
www.ukmedical.com

UKMEDICAL[®]
It's Interventional.