

UKMEDICAL®

*It's Interventional.*

*It's about expanding  
embolisation  
solutions /*



**IMPEDE® & IMPEDE-FX**

[hello@ukmedical.com](mailto:hello@ukmedical.com)  
[ukmedical.com](http://ukmedical.com)




# *It's about peripheral vascular embolisation /*

*The IMPEDE® Embolisation Plug is indicated to obstruct or reduce the rate of blood flow in the peripheral vasculature.*

*Comprised of a self-expanding porous shape memory polymer, IMPEDE® is crimped for delivery and undergoes a heat driven shape change once deployed in a vessel. Rapid thrombosis occurs creating an effective, stable occlusion.*

## **Product benefits:**


- *High embolic material volume*
- *Efficient, predictable space filling*
- *Easily navigates tortuosity*
- *Minimal artefact*
- *Accelerated healing*
- *Simple to use*
- *Accurate placement*
- *Stable occlusion*

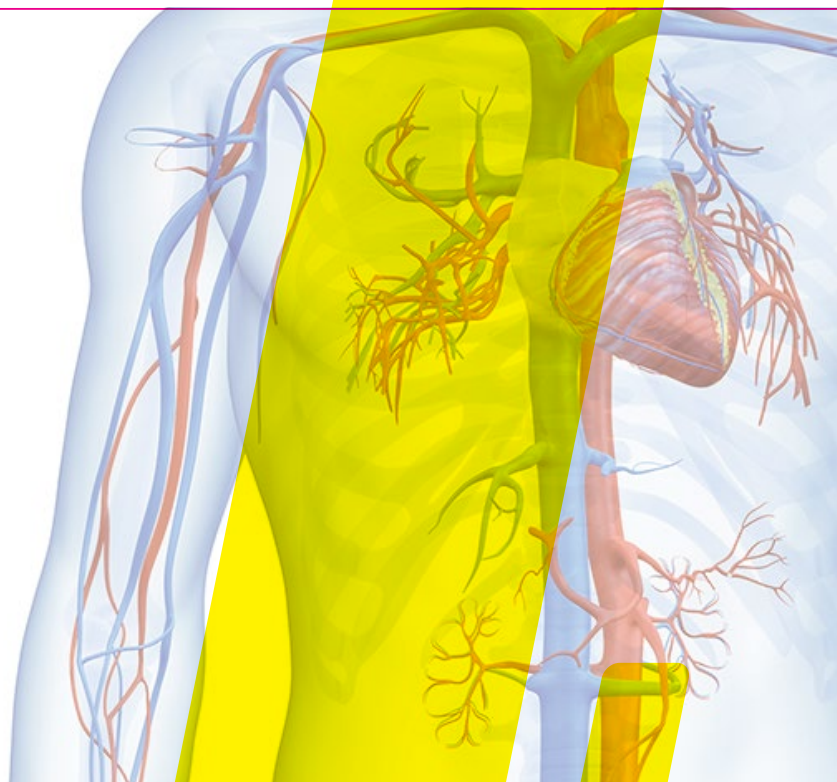

## **Additional information?**

If you would like further assistance with using our Shape Memory range, or need advice on how to order, please call us on: **0114 268 8880**


**Product features:**

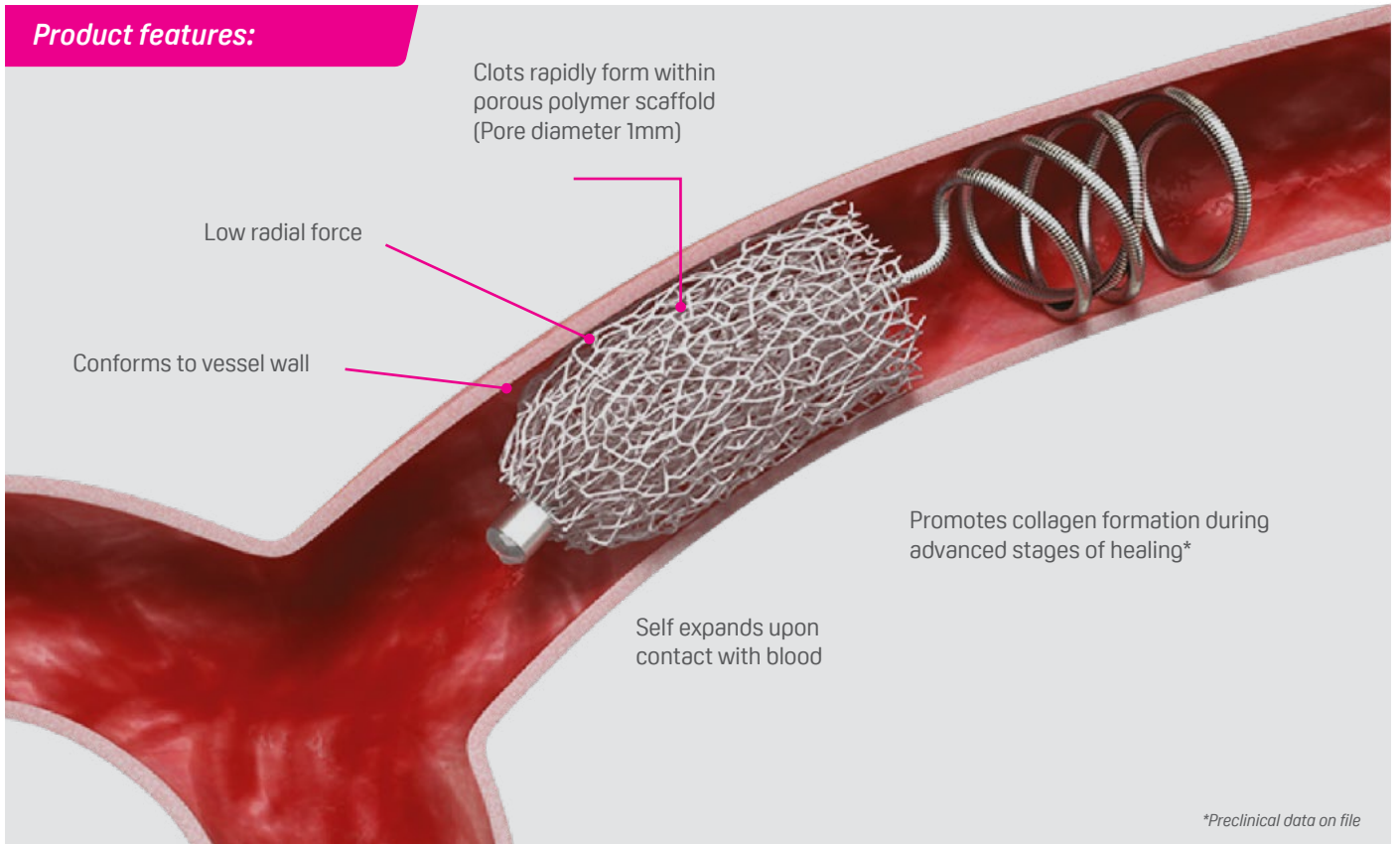


**Delivery system:**



**IMPEDE®**

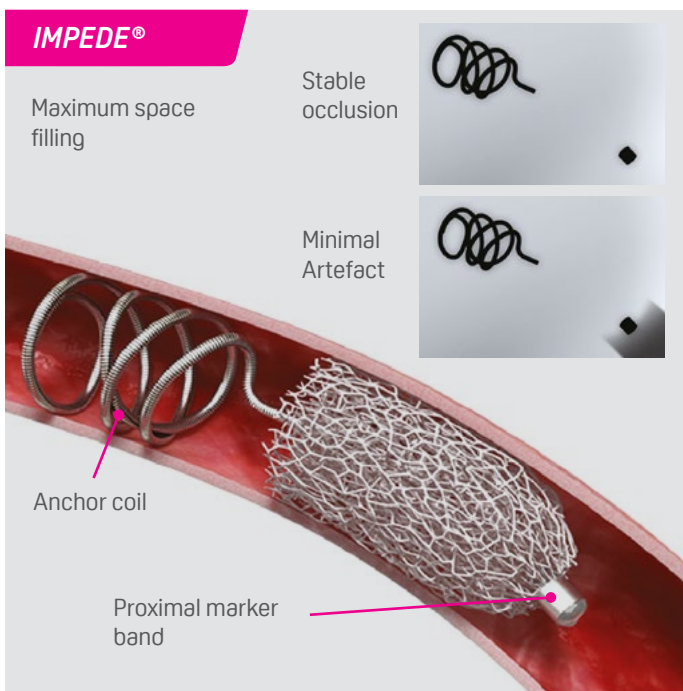

**IMPEDE-FX:**

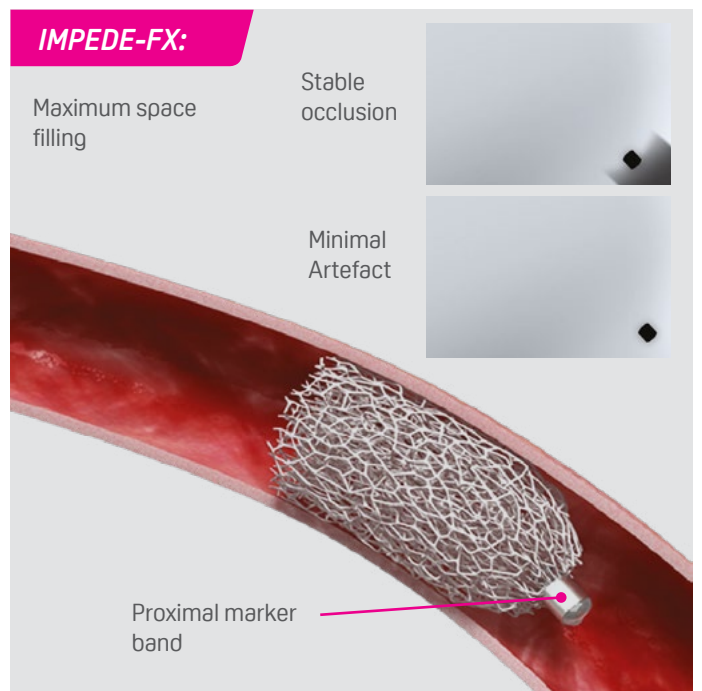

## Ordering Information

IMPEDE®

IMPEDE-FX



Description	Product code	Vessel dia. mm	Proximal marker dia. (inch / mm)	Expanded SMP plug dia. (mm)	SMP plug length	Anchor coil dia. (mm)	Catheter ID* min/ max (mm)	Catheter ID* min/ max (inches)
IMPEDE®	IMP-05	2-5	.032" 0.81mm	6	10	7	0.97 / 1.45	.038 / .057
	IMP-07	4-7	.046" 1.17mm	8	10	9	1.40 / 1.80	.055 / .071
	IMP-10	6-10	.065" 1.65mm	12	15	13	1.78 / 2.29	.070 / .090
IMPEDE-FX	IMP-FX-06	-	.032" 0.81mm	6	10	-	0.97 / 1.45	.038 / .057
	IMP-FX-08	-	.046" 1.17mm	8	10	-	1.40 / 1.80	.055 / .071
	IMP-FX-12	-	.065" 1.65mm	12	15	-	1.78 / 2.29	.070 / .090

\*Catheter recommendations based on .035" guidewire.



Distributed by:

UK Medical  
Albreda House  
Lydgate Lane  
Sheffield, S10 5FH

T: +44 (0) 114 268 8880  
F: +44 (0) 114 268 8881  
sales@ukmedical.com  
www.ukmedical.com

@UK\_Medical

